






**UNDERBODY**  
SERIES

# Bair Hugger® Therapy Lithotomy Underbody Blanket

Model 585



Only from  
Arizant Healthcare Inc.

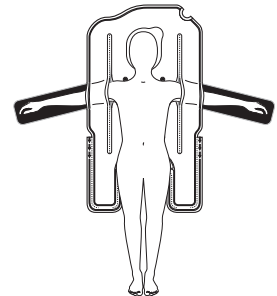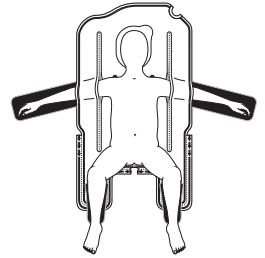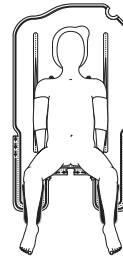

The Lithotomy Underbody Blanket is positioned on the procedure table before the patient arrives to the room.

Designed to deliver full, unrestricted patient access and the flexibility for procedures involving the lower extremities and abdominal, peritoneal and pelvic cavities.



# Lithotomy Underbody Blanket Features Model 585


**UNDERBODY**  
SERIES



- Underbody design provides full unrestricted access to the patient
- Fluid outlets minimize pooling of fluids on the surface of the blanket
- Pass-through slit design allows for flexible patient positioning and accommodates the use of a drawsheet
- Underbody design accommodates lateral, supine and lithotomy patient positions
- Generously-sized peritoneal cut-out at the base of the blanket provides clinician unobstructed access to patient
- Integrated tie strips can be loosely secured around the stirrups/leg supports
- The blanket draws up next to patients when inflated to surround them in warmth
- Tape and tuck flaps secure the blanket to the OR table



## Product Offerings

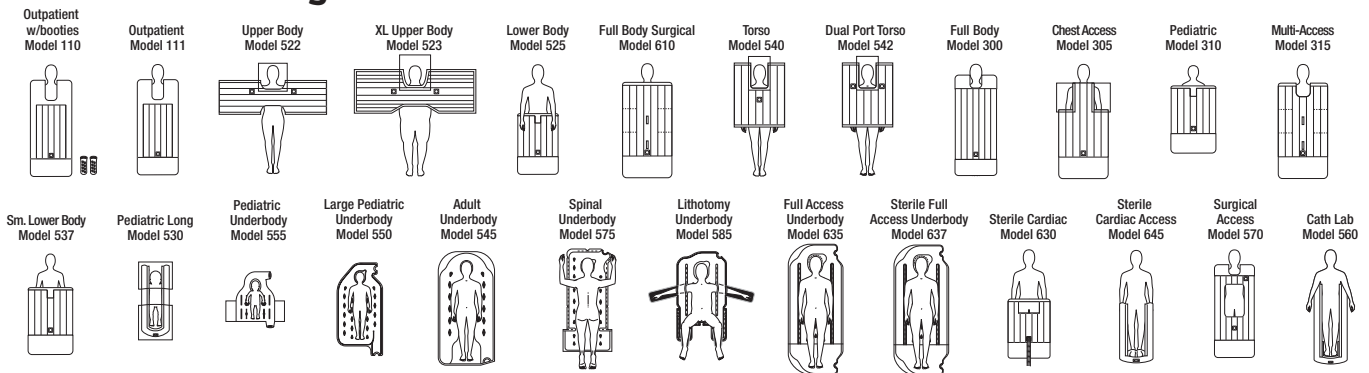

©2008 Arizant Healthcare Inc. All rights reserved.  
 10393 West 70th Street, Eden Prairie, MN 55344 USA  
 tel: 952-947-1200 | 800-733-7775  
 fax: 952-947-1400 | 800-775-0002  
 web: arizanthealthcare.com | bairhugger.com

Bair Hugger, Arizant, Arizant Healthcare, and the Bair Hugger and Arizant logos are trademarks of Arizant Healthcare Inc., registered or pending in the U.S. Patent & Trademark Office and in other countries.

601907D 11/08

### Specifications

**Lithotomy Underbody Blanket - Model 585**  
 Size: 64" x 36" (160 x 91 cm)  
 Weight: 5.4 oz (153 g)  
 Drape: 24" x 24" (61 x 61 cm)

### Order Information

Part Number: 58501  
 Units/Case: 10

To order, please contact your local Arizant Healthcare sales representative or call Customer Service at 1-800-733-7775.