


RANGER® SYSTEM PEDIATRIC FLUID WARMING

Ranger Pediatric/Neonate Blood and Fluid Warming Disposable Sets

The unique design of the pediatric/neonate disposable sets' heat exchanger allows the Ranger system to fully address pediatric fluid warming needs. These DEHP-free sets are designed for low flow, maintenance fluid delivery and aliquot fluid administration where measured fluid volumes are infused via syringe. In addition, the priming volume of these sets is 20mL. Blood and IV fluids can be warmed up to six liters per hour.


The Ranger system's pediatric/neonate sets complement the Bair Hugger® product line's multiple styles of pediatric warming blankets to help healthcare facilities effectively manage their pediatric warming needs.

Arizant Healthcare is committed to designing and manufacturing safe and effective patient warming products for all ages and for a variety of procedures.

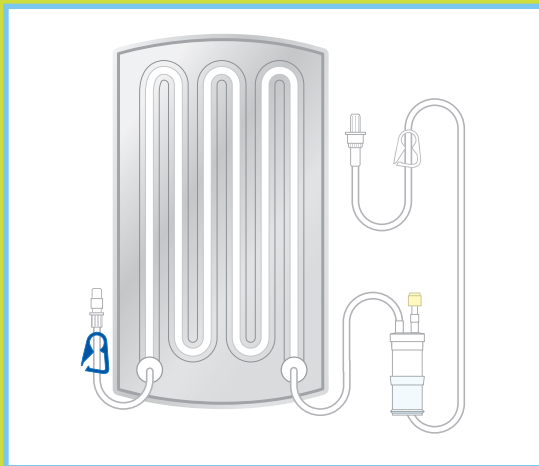

Pediatric/Neonate Disposable Set
Model 24400

Flow rates: KVO - 6,000 mL/hr
Priming volume: 20 mL
Components: Heat exchanger, bubble trap with air aspiration port, 24 in. patient line
Pressure infusion: Up to 300 mmHg
Material: Latex-free, DEHP-free
Case quantity: 10 disposable sets/case
Sterile: EtO
Use: Single use only

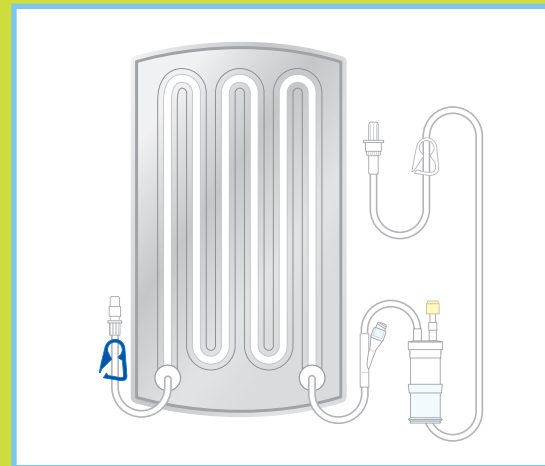

Pediatric/Neonate Disposable Set with Fluid Aspiration Port / Model 24450

Flow rates: KVO - 6,000 mL/hr
Priming volume: 20 mL
Components: Heat exchanger, bubble trap with air aspiration port, fluid aspiration port, 24 in. patient line
Pressure infusion: Up to 300 mmHg
Material: Latex-free, DEHP-free
Case quantity: 10 disposable sets/case
Sterile: EtO
Use: Single use only



Blood and Fluid Warming Systems

www.rangerfluidwarming.com

©2010 Arizant Healthcare Inc. All rights reserved. 10393 West 70th Street, Eden Prairie, MN 55344
 USA tel: 952-947-1200/800-733-7775 fax: 952-947-1400/800-775-0002. Ranger, Bair Hugger
 and the Bair Hugger logo are trademarks of Arizant Healthcare Inc., registered in the U.S. Patent and
 Trademark office. 601479C 8/10